



THE UNIVERSITY *of* EDINBURGH

Edinburgh Research Explorer

The Valgent4 protocol

Citation for published version:

Bonde, J, Ejegod, DM, Cuschieri, K, Dillner, J, Heideman, DAM, Quint, W, Pavon Ribas, MA, Padalko, E, Christiansen, IK, Xu, L & Arbyn, M 2018, 'The Valgent4 protocol: Robust analytical and clinical validation of 11 HPV assays with genotyping on cervical samples collected in SurePath medium', *Journal of Clinical Virology*, vol. 108, pp. 64-71. <https://doi.org/10.1016/j.jcv.2018.09.012>

Digital Object Identifier (DOI):

[10.1016/j.jcv.2018.09.012](https://doi.org/10.1016/j.jcv.2018.09.012)

Link:

[Link to publication record in Edinburgh Research Explorer](#)

Document Version:

Publisher's PDF, also known as Version of record

Published In:

Journal of Clinical Virology

General rights

Copyright for the publications made accessible via the Edinburgh Research Explorer is retained by the author(s) and / or other copyright owners and it is a condition of accessing these publications that users recognise and abide by the legal requirements associated with these rights.

Take down policy

The University of Edinburgh has made every reasonable effort to ensure that Edinburgh Research Explorer content complies with UK legislation. If you believe that the public display of this file breaches copyright please contact openaccess@ed.ac.uk providing details, and we will remove access to the work immediately and investigate your claim.





The Valgent4 protocol: Robust analytical and clinical validation of 11 HPV assays with genotyping on cervical samples collected in SurePath medium

Jesper Bonde^{a,*}, Ditte Møller Ejegod^a, Kate Cuschieri^b, Joakim Dillner^c, Daniëlle A.M. Heideman^d, Wim Quint^e, Miguel Angel Pavon Ribas^f, Elizaveta Padalko^g, Irene Kraus Christiansen^h, Lan Xuⁱ, Marc Arbynⁱ

^a Molecular Pathology Laboratory, Department of Pathology, Copenhagen University Hospital, Kettegårds Alle 30, 2650 Hvidovre, Denmark

^b Scottish HPV Reference Laboratory, Royal Infirmary of Edinburgh, EH16 4SA Edinburgh, Scotland, United Kingdom

^c International HPV Reference Center, Karolinska University Hospital, Alfred Nobels Allé 8, Huddinge, Stockholm, Sweden

^d Cancer Center Amsterdam, Department of Pathology, VU University Medical Center, De Boelelaan 1118, 1081 HV Amsterdam, The Netherlands

^e DLL Diagnostic Laboratory, Visseringlaan 25, 2288 ER Rijswijk, The Netherlands

^f Infection and Cancer Laboratory, Cancer Epidemiology Research Program and CIBER-ONC, Granvia de L'Hospitalet 199-203, Institut Català d'Oncologia, Barcelona, Spain

^g Department of Laboratory Medicine, Ghent University Hospital, C. Heymanslaan 10, B-9000 Ghent, Belgium

^h Norwegian HPV Reference Laboratory, Department of Microbiology and Infection Control, Akershus University Hospital, 1478 Lørenskog, Norway

ⁱ Sciensano, Belgian Cancer Centre/Unit of Cancer Epidemiology Juliette Wytsmanstreet 14, B1050 Brussels, Belgium

ARTICLE INFO

Keywords:

HPV genotyping assays
VALGENT
Cervical cancer screening
Clinical validation

ABSTRACT

Background: The VALidation of HPV GENotyping Tests (VALGENT) is an international initiative designed to validate HPV assays with genotyping capability. The VALGENT4 protocol differs from previous VALGENT installments as the sample collection medium is SurePath, and exclusively includes samples from women ≥ 30 years of age which is concordant with the majority of HPV primary screening guidelines. Here we present the protocol for the fourth installment of the VALGENT framework.

Objectives: In VALGENT4 11 HPV assays will be evaluated using two comparator assays based on PCR with the GP5+/6+ primers.

Study design: Overall, the VALGENT4 panel consists of 1,297 routine samples comprised of 998 unselected, consecutive samples, of which 51 samples had abnormal cytology with 13 women diagnosed with \geq CIN2, and 299 consecutive samples enriched for \geq ASCUS cytology (100 ASCUS, 100 LSIL, 99 HSIL) with 106 \geq CIN2 upon follow up. Manipulated and DNA extracted panel samples were characterized with respect to human beta globin (HBB) and overall DNA content and composition to quality assess the panel prior to distribution to the collaborating sites.

Result: The relative cellularity (mean CT value of HBB from the Onclarity assay) on the 1,297 LBC samples (CT=24.8) was compared with 293 un-manipulated routine cytology screening samples (CT=23.8). Furthermore, the DNA extracted panel samples was characterized using the Exome iPLEX pro assay, which reports amplifiable copies on individual samples as well as copies of five different base pair lengths. Here the data showed a slightly lower number of amplifiable DNA copies (ratio: 0.7, $p = < 0.01$) in the VALGENT4 panel samples compared to routine extracted cervical DNA samples

Conclusion: The present manuscript details the manipulation, processing and quality assessment of samples used in VALGENT-4. This methodological document may be of value for future international projects of HPV test validation.

1. Background

Practically all cervical cancers are derived from an Human papillomavirus (HPV) infection [1–3] and the causal relation is firmly

established between the development of cervical cancer and at least 12HPV genotypes (HPV 16, 18, 31, 33, 35, 39, 45, 51, 52, 56, 58, and HPV59, all group 1 carcinogens as classified by IARC) [1,2]. Evidence on the clinical value of HPV testing for triage of borderline cervical

* Corresponding author.

E-mail address: jesper.hansen.bonde@regionh.dk (J. Bonde).

<https://doi.org/10.1016/j.jcv.2018.09.012>

Received 20 June 2018; Received in revised form 7 September 2018; Accepted 14 September 2018

1386-6532/ © 2018 The Authors. Published by Elsevier B.V. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

cytology and test of cure in women who have been treated for pre-invasive cervical lesions is today widely accepted [4–6]. Additionally, randomized controlled trials have provided evidence for the use of HPV-based screening over cytology as a superior and highly reproducible screening modality to prevent both invasive squamous and adeno-carcinomas [5,7]. As a result, primary HPV-based cervical screening is being implemented in several countries globally, substituting conventional and liquid based cytology (LBC) as the standard of care. Together, this has created a rapidly expanding market for molecular HPV tests, and more than 200 different HPV assays are now commercially available [8]. As primary HPV based screening becomes the standard of care, additional focus on evidence of clinical performance of HPV assays for use in cervical cancer screening programs is required.

HPV assays, applicable to cervical screening and disease management contexts, have undergone a rapid technical and scientific evolution over the last decade. The first generation of clinical HPV assays were developed solely to detect oncogenic HPV genotypes using DNA PCR or hybridization techniques and mainly reported the test outcome as either HPV positive or negative, with no individual HPV genotype information available. Newer generations of commercially available HPV assays depend upon a wide variety of DNA and RNA detection techniques and allow for individual reporting of HPV genotypes. Current HPV genotyping assays can basically be divided into three categories (1) Assays with limited genotyping that report separate identification of HPV16 and HPV18 or HPV18/HPV45, combined with pooled detection of the remainder of the oncogenic types, (2) Assays with extended genotyping that report separate identification of ≥ 5 genotypes combined with one or more bulk detections of the remainder, and (3) Assays with full genotyping, reporting individual identification of all carcinogenic HPV genotypes [9].

Studies have assigned different carcinogenic potential, or risk, to individual HPV genotypes, with HPV16 being the most potent, in particular for squamous carcinomas of the cervix [10–15], followed by HPV18 and to some extent HPV45 associated with adenocarcinomas of the cervix [12,16–18]. Together, up to 70% of all cervical cancers are caused by HPV16 and HPV18. An additional five HPV genotypes (HPV31, 33, 45, 52, & 58) add around 19% to the burden of cervical cancer incidence [16]. Recent study data from long term follow up of a Danish cervical screening cohort, have shown that persistent infections with HPV16, HPV18, HPV31 and HPV33 are associated with an increased risk for development of cervical intraepithelial neoplasia grade 3 or worse (\geq CIN3) compared to the remaining oncogenic genotypes [14,19]. The six carcinogenic types HPV35, 39, 51, 56, 59, 68 in contrast contribute 8–9% of the cumulative proportion of HPV-induced cervical cancer [16,20]. Other HPV genotypes incl. HPV66 and 68 are only rarely involved in cervical carcinogenesis, but both are included in all oncogenic HPV assays [16]. Finally, genotype distribution among highly vaccinated cohorts of women have drastically changed affecting the performance of both cytology screening and HPV genotyping [21].

Taken together, knowledge of the genotype(s) in a screening sample may add information on the individual woman's risk of developing disease [12]. Consequently, HPV technology is progressing from simple plus/minus outcomes towards revealing more detailed information on the HPV infection in question. This fact is reflected in a growing number of national screening guidelines where HPV16/HPV18 genotyping is integrated in screening algorithms [5,22,23].

Fulfillment of the “Guidelines for human papillomavirus DNA test requirements for primary cervical cancer screening in women of 30 years and older” [24], is widely accepted as the methodology by which to validate HPV assays for screening use [24,25] and several limited or extended genotyping HPV assays have fulfilled the international guidelines [25–32]. The international validations are based upon the overall clinical performance but not with respect to the clinical performance of detection of individual genotypes *per se*. Assay specific detection of individual genotypes has instead been addressed in studies of defined

plasmid based panels with known copy numbers [33,34]. The VALi-dation of HPV GENotyping Tests (VALGENT) initiative was established to create a framework for comparison and validation of HPV genotyping assays using validated comparator assays, similar to the international guidelines [24] and using a relevant sample population with sufficient disease to confirm clinical performance while including primary screening samples [9]. So far, three installments of VALGENT protocols have occurred, where a broad variety of HPV assays with both extended and full genotyping capability have been evaluated [9].

2. Objectives

The fourth installment of the VALGENT framework, the VALGENT4 study, aims to provide high quality comparative data on clinical performance of HPV assays on cervical screening samples. The VALGENT4 protocol differs from previous VALGENT installments as the sample collection medium is BD SurePath™ (Becton, Dickinson and Company), targets DNA assays only, and exclusively includes only samples from women ≥ 30 years of age which is concordant with the majority of HPV primary screening guidelines. In addition, a high number of assays previously untested on this LBC collection media will be evaluated. Moreover, the VALGENT4 protocol is unique in the VALGENT iterations by proposing a new standard for describing the quality of the sample panel included.

3. Study design

3.1. Study description

Sample collection, registration, processing and aliquoting was done at the parent site at the Molecular Pathology Laboratory, Dept. Pathology, Copenhagen University Hospital, Hvidovre. The VALGENT4 samples will be provided to participating scientific partners anonymized in concordance with Danish Data Protection Agency regulations. Aliquot panels were shipped refrigerated from parent laboratory to scientific partners.

The VALGENT4 panel was generated using fresh SurePath collected screening samples from the Danish cervical cancer screening program which services a well-screened population with a high background risk of cervical cancer [35]. Moreover, the parent laboratory delivered the panel in two defined versions; pre-extracted and quality described DNA for assays requiring DNA as input material (time between collection and DNA extraction: Mean 28 days, min: 2 days, max: 70 days), and original LBC material (time between collection and aliquoting for panel: Mean 16 days, min: 2 days, and max: 67 days) to those laboratories where assays with full, integrated CE-IVD marked work flows were interrogated. The option to have a panel consisting of extracted DNA limits any variability of various 3rd party DNA extraction platforms on HPV assay performance. Finally, a novel element entailed characterization of all included panel samples using a MALDI-TOF assay which measured the available DNA as well as the relative level of DNA fragmentation in the individual sample.

In total, eleven different HPV genotyping assays from 8 different manufacturers will be evaluated in the VALGENT4 study (Table 1), using GP5+/6+ PCR-EIA with genotyping as comparator in line with previous VALGENT installments.

3.2. Sample collection

The VALGENT4 study was embedded into the routine cytology-based operations of the Department of Pathology, Copenhagen University Hospital, Hvidovre, Denmark. This pathology department receives virtually all cervical cytology screening samples from women residing in the Danish Capital Region, covering almost a third of the entire Danish population. Women are invited for screening every three years at ages 23–49 and every fifth year at ages 50–59. An HPV based

Table 1

HPV genotyping assays evaluated, concurrent material required and scientific partners under the VALGENT4 study protocol.

Valgent-4 included Assays	Aliquot	Amplicon length	Scientific partner
BD Onclarity HPV Assay	Original Material	79–137 bp.	Hvidovre Hospital, Dept. Pathology, Hvidovre, Denmark
Genomica CLART HPV4 assay	DNA	465 bp.	
Agena HPV MassArray assay	DNA	90–122 bp.	
Roche cobas 4800 HPV Test	Original Material	~ 200 bp	Norwegian HPV Reference Lab, Akershus University Hospital Norway
Fujirebio INNO-LiPA Genotyping Extra II test	DNA	65 bp	Ghent University Hospital, Ghent, Belgium
SeeGene Anyplex HPV28 detection test	DNA or Original	~ 150 bp	Infection and Cancer Laboratory. Cancer Epidemiology Research Program, Institut
SeeGene Anyplex II HPV test	Material	~ 150 bp	Català d'Oncologia Barcelona, Spain
Self-screen HPV-Risk assay	DNA	~ 150 bp	Cancer Center Amsterdam, VU University Medical Center, Amsterdam, The Netherlands
GeneFirst HPV-MPA Genotyping Test	DNA	150 bp.	Scottish HPV Reference Lab,
Liferiver Harmony test	DNA	100–200 bp.	Royal Infirmary of Edinburgh
Liferiver Venus tes		100–200 bp.	Scotland
Comparator assays			
GP5+/6+ EIA Luminex	DNA		DDL Diagnostic Laboratory, Rijswijk, The Netherlands
GP5+/6+ PCR EIA kit HPV GP HR	DNA		International HPV Reference Center, Karolinska University Hospital, Stockholm, Sweden

“check out test” is offered to women 60–64 years of age prior to exiting the screening program at age 65. The Danish cervical cancer screening program is free of charge at point of delivery. Samples are collected predominately by general practitioners or in some cases by gynecologists. The Combi-brush (Rovers, Oss, Netherlands) is used for collection of cell material and all cytology samples are collected in SurePath medium.

Samples from women aged 30–59 years were included in VALGENT4 study.

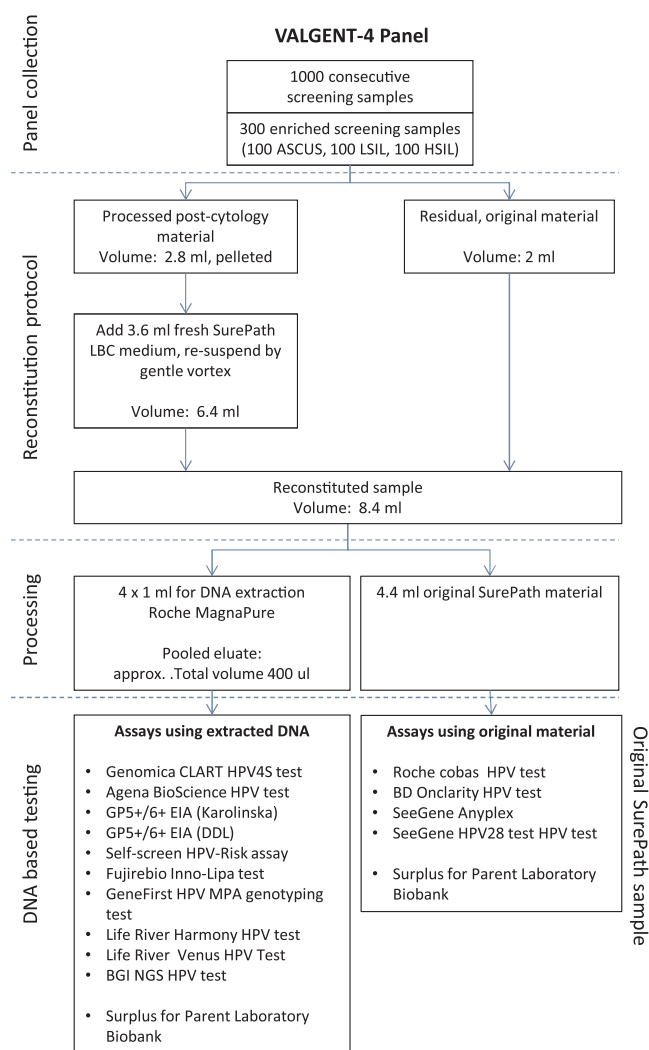
The VALGENT4 panel consists of 1297 routine samples collected in two groups as previously described by the VALGENT framework [9]:

- 998 unselected, consecutive samples, of which 51 samples had abnormal cytology (\geq atypical squamous cells of undetermined origin, (ASCUS)). By subsequent histological follow-up, 13 women were diagnosed with cervical intraepithelial neoplasia 2 or more (\geq CIN2)
- 299 consecutive routine samples enriched for \geq ASCUS cytology (100 ASCUS, 100 low-grade squamous intraepithelial lesions (LSIL), 99 high-grade squamous intraepithelial lesions (HSIL)). From this population, 106 women were presented with \geq CIN2 upon follow up.

Three samples were discarded after collection due to lack of general consent after mandatory cross-check in *Danish human biological material in health research projects register* (Vævsanvendelsesregisteret).

3.3. Sample processing

SurePath samples nominally contain 10 ml of liquid upon arrival in the laboratory. Currently, the Department of Pathology processes cervical cytology samples using the BD Totalys instrumentation. Here, the cytology processing uses 8 ml of the original sample, which is column purified for debris and mucus, spun down, pelleted, and re-suspended in 1 ml BD Density Reagent. Of this 1 ml medium, 200 μ l is used for the cytology slide procedure. In routine, the remaining 800 μ l is automatically added the 2 ml of fresh SurePath medium, resulting in a residual volume of 2.8 ml. By calculation, the final volume of the pelleted material in the 2.8 ml holds a 2.3 times higher concentration of cellular material than the original material. In order to approximate the cellularity of the original 10 ml sample, we developed a reconstitution protocol utilizing post-cytology pelleted material mixed with residual original material. The post-cytology pelleted sample was reconstituted by adding 3.6 ml of fresh SurePath medium, resulting in a volume of 6.4 ml (Fig. 1). In addition, the residual, original 2 ml surplus material from the cytology procedure was added to the 6.4 ml of reconstituted

**Fig. 1.** Flowchart.

sample. Consequently, the nominal volume of the reconstituted sample available for testing in the VALGENT4 protocol was 8.4 ml.

The reconstituted samples were split into aliquots for DNA extraction and original sample material upon processing of the samples (see processing time above).

3.4. DNA extraction

DNA was extracted under strictly controlled conditions using standard operating procedure clinically routine DNA extraction protocol at the parent laboratory based upon the MagNA Pure96 system (Roche Diagnostics, Pleasanton, CA). In short, 4 × 1 ml aliquots of reconstituted material were transferred to four Eppendorf tubes, spun down (14,000 rpm, 5 min), the supernatant removed, and the pellet re-suspended in a mix of 180 µl phosphate buffered saline (10x conc. pH 7.4, Pharmacy product) and 20 µl Proteinase K (Recombinant, PCR grade, Roche Diagnostics). The samples were vortexed and incubated one hour at 56 °C for Proteinase K digestion. Subsequently, the tubes were incubated for one hour at 90 °C to reverse SurePath formaldehyde induced co-valent cross linking. The entire volume was transferred to the MagNA Pure 96 system and extracted using the MagNA Pure LC Total Nucleic Acid Isolation Kit (Roche Diagnostics). The four DNA eluates were pooled into one tube to a total volume of 400 µl.

3.5. Cytology

Cytology was reported according to the Bethesda 2001 criteria, using computer assisted reading by BD Focal Point and Slide Wizard systems. ASCUS cases were triaged (reflex testing) routinely with HPV testing (BD Onclarity). All abnormal cytology findings were adjudicated as per standard operating procedures of the parent laboratory by a senior pathologist. Women with LSIL were recommended for repeated cytology testing after 6 months. Women with inadequate cytology were recommended for repeated sampling and testing within 3 months. Women with normal cytology were returned to routine screening according to age-specific intervals. The complete screening history of the women included in VALGENT4 was retrieved from the Danish Pathology Data Bank (PatoBank). All clinical follow up was managed according to Danish guidelines, and the outcome of the VALGENT4 HPV testing did not affect clinical follow up recommendations. All HPV testing performed in the framework of VALGENT4 was blinded to cytology and clinical follow-up.

3.6. Histology

All histology available for disease ascertainment was derived from the clinical follow-up of women whose samples were included in the VALGENT4 cohort. The following screening sample outcomes elicit a recommendation for follow-up with colposcopy and biopsy-taking under the current Danish guidelines:

- Women, ≥30 years with ASCUS and HPV-positivity upon triage testing.
- Women with HSIL, atypical squamous cells – cannot exclude HSIL (ASC-H), atypical glandular cells (AGS) or cytological indications of carcinoma.
- Women with continued ASCUS and LSIL cytological diagnosis, as evaluated by the present sample and screening history.

Danish Gynecology Guidelines recommend biopsy taking on all aceto-white lesions observed, or a random four quadrants biopsy set where no lesions are visible upon colposcopy. The histological follow-up for all women included in the VALGENT4 study was retrieved from the PatoBank.

3.7. Data source for clinical outcomes

The PatoBank follow-up period was on average 19 months (min. 18 months; max 20 months) from collection of the VALGENT4 samples. The registration of cervical cytology and histology in the PatoBank has been complete on a national level since 2008.

3.8. Participating scientific partner panel testing

All HPV testing within VALGENT4 was blinded to cytology and histology outcomes from the parent laboratory as well as the HPV test results from other scientific partners. All test results were reported to the Scientific Institute of Public Health, Brussels, where the data are compiled, and performance calculations are conducted [9]. All HPV testing was performed according to manufacturer's specifications. Nine different participating scientific partners (including the parent site) are taking part in the VALGENT4 study (Table 1).

3.9. Pre-analytical characterization of LBC material and DNA aliquot panel samples

3.9.1. Pre-analytical characterization of cellularity of reconstituted original LBC material

All 1297 panel samples were tested with the BD Onclarity HPV assay (Onclarity) on the automated Viper Lt platform under the VALGENT4 protocol. The Onclarity assay harbors an internal Human Beta globin (HBB) control for assay performance and sample adequacy. To evaluate the resulting cellularity of the resuspension protocol, mean HBB Ct values from the 1297 VALGENT4 samples were retrieved from the Viper Lt and used as a proxy marker for cellularity. As comparator, mean HBB Ct values of unselected, unannotated routine cervical screening SurePath samples were used as a reference group for cellularity (designated "Control-1", N = 293). The Onclarity assay on the Viper Lt Platform workflow has previously been described in detail [27,36]. In summary, 0.5 ml original SurePath material is transferred to a BD CBD medium tube prior to heat treatment for 30 min at 120 °C on the BD pre-warm station. The pre-warmed samples are subsequently transferred to the fully automated Viper Lt platform and tested with the Onclarity assay according to manufacturer's recommendations. Mean HBB values between the two groups were compared using one-way ANOVA testing (IBM SPSS ver. 22).

3.9.2. Amplifiable DNA copies and DNA fragmentation evaluation of DNA extracts

The iPLEX PRO ExomeQC assay (Agena Bioscience, Hamburg) is a quantitative assay which reports average amplifiable human DNA copies in a sample, in addition to reporting the relative number of amplifiable fragments of 100, 200, 300, 400, and 500 base pairs (bp) lengths. Here, the assay was used as a quality control of the DNA extracts derived from the reconstituted VALGENT4 panel. As comparator, DNA from unselected, unannotated, and un-manipulated SurePath cervical screening samples collected from the routine service were used as reference (designated "Control-2", N = 184).

The ExomeQC panel harbors 44 Single Nucleotide Polymorphisms (SNPs), three markers for gender identification and five copy-number markers in a single multiplexed assay. The assay uses the Matrix Assisted Laser Desorption/Ionization-Time of Flight Mass Spectrometry (MALDI-TOF MS) technology. An initial multiplex PCR amplification was set up with 2 µl DNA (duration: 2.5 h), followed by a Shrimp Alkaline Phosphatase (SAP) reaction (dephosphorylize excess Nucleotides, duration: 50 min). This was followed by an iPLEX Pro single base extension PCR reaction in which a mix of oligonucleotide extension primers designed to anneal to the amplified DNA fragments, was added together with an extension enzyme and mass-modified dideoxynucleotide terminators (duration 3.5 h). The extension products were subsequently desalted with Clean Resin prior to being loaded into the MassArray Dx Nanodispenser RS1000 (Agena Bioscience, Hamburg), which transfers the analyte to a spectroCHIP. Here the sample crystallize with the matrix on the chip, which was then analyzed on the MassARRAY Dx Analyzer 4 (Agena Bioscience, Hamburg). The crystals were irradiated by a laser, inducing desorption and ionization. The MA4 accelerates the samples to a detector that differentiates genetic variants by molecular mass.

Table 2
Overall characteristics of the VALGENT4 study population.

Characteristics	Screening population		Enriched population ^a		Total population	
	N	%	N	%	N	%
Age						
Mean Age in years	42.8 (36.0–49.0)		40.4 (34.0–45.8)		42.2 (36.0–48.0)	
30–39	383	38.4%	148	49.8%	531	41.0%
40–49	408	40.9%	111	37.4%	519	40.1%
50–59	207	20.7%	38	12.8%	245	18.9%
30–59	998	100%	297	100%	1,295	100%
Cytology						
NILM	947	94.9%	0	0%	947	73.0%
ASCUS	6	0.6%	100	33.3%	106	8.2%
LSIL	21	2.1%	100	33.3%	121	9.3%
HSIL. AGC. ASC-H	24	2.4%	99	33.3%	123	9.5%
Total	998	100%	299	100%	1297	100%
Histology						
No histology	911	91.3%	55	18.4%	966	74.6%
Normal	57	5.7%	82	27.4%	139	10.7%
CIN1	17	1.7%	54	18.1%	71	5.5%
CIN2	5	0.5%	34	11.4%	39	3.0%
CIN3	7	0.7%	64	21.4%	71	5.5%
Cancer	1	0.1%	8	2.7%	9	0.7%
All histologies	87	8.7%	242	80.9%	329	25.4%
< CIN2	74	7.4%	136	45.5%	210	16.2%
≥ CIN2	13	1.3%	106	35.5%	119	9.2%

^a Two women enrolled in Valgent4 each had two samples included; therefore therefor 299 samples from 297 women are included in the enriched cohort.

3.9.3. Statistical analysis

The Onclarity assay has a three well set-up with nine HPV genotype read-outs with internal HBB control included in each well. The Ct-value of the HBB was calculated as an average of the three HBB Ct values in each of the three wells for every panel and control sample. The Onclarity cut-off for all channels is Ct 34.2. For the ExomeQC analysis, inadequate samples (insufficient number of SNPs detected) and samples with outlier values above 10,000,000 units were excluded in the resulting analysis. Mean Ct, mean available copy numbers, 95% confidence interval, standard deviations as well as the one-way ANOVA test used to calculate the statistical difference in cellularity and DNA integrity parameters between the VALGENT4 panel and Control-1 and 2 samples were calculated using IBM SPSS statistics 22.

4. Results

The overall characteristics of the VALGENT4 study population are summarized in Table 2. The mean age was 42.8 years for all included women, and 40.4 and 42.2 years for the unselected, consecutive and enriched subpopulations, respectively. In the 998 consecutive, unselected panel samples, 5.1% (51 samples) had abnormal cytology and as of March 2018, 1.3% (13 samples) had ≥ CIN2 follow-up histology registered in the Danish National Pathology Registry (PatoBank). The enriched population of 100 ASCUS, 100 LSIL and 99 HSIL cytology samples resulted in 106 ≥ CIN2 histology samples as per register update March 2018. In total, 119 ≥ CIN2 cases were registered at data retrieval ≥ 20 month after baseline.

Comparing the relative cellularity (mean Ct values of HBB control from the Onclarity assay) on the 1,297 LBC panel samples with data from 293 separately collected, un-manipulated routine cytology screening samples (Denominated “Control-1” samples, Table 3A). The mean Ct value of the VALGENT4 samples (Ct = 24.8) was found to be on average 1.0 Ct above the mean Ct value of the Control-1 samples (Ct = 23.8), indicating slightly less cellularity in the panel samples compared to native samples.

For DNA extracted panel samples, characterization was done using the Exome iPLEX pro assay, which reports amplifiable copies on individual samples as well as amplifiable copies at five different base pair

Table 3A

Quality assurance of the VALGENT4 LBC samples: Mean HBB Ct value as approx. marker for cellularity.

Group	N	Mean HBB Ct	Ratio ^a	St.d	95% Confidence Incidence		Min	Max	P-value ^b
					Lower	Upper			
Valgent4	1297	24.8	1.04	1.5	24.7	24.9	20.8	32.3	0.00
Control1	293	23.8	1.00	1.5	23.6	24.0	20.5	32.0	

^a The Control1 population was used as a reference.

^b p-value is calculated using the one-way ANOVA test.

lengths (Table 3B). Here, the data showed statistically difference in the number of amplifiable DNA copies comparing the VALGENT4 panel to the Control-2 population (ratio of 0.7, p-value < 0.01). When looking at the relative proportion of different base pair length fragments, the statistical difference was only evident at the longest 400 and 500 bp fragments (Table 3B).

5. Discussion

The VALGENT framework is an international cooperation aimed at evaluating HPV genotyping tests for clinical use in cervical cancer screening. Here we present the study protocol for the fourth installment of the VALGENT framework, presenting data that provides insight into the underlying quality of the 1,297 VALGENT4 samples. One of the novel elements in VALGENT4 is the resuspension protocol employed which attempts to emulate the original cellularity of un-manipulated screening samples (Fig. 1).

For research use, most cross-sectional studies where HPV tests are applied rely on stored, residual clinical cytology screening material. In contrast, most HPV tests are intended to run on recently collected and un-manipulated original sample material. In addition, for evaluations where multiple assays are applied to samples, the volume of sample may be a rate-limiter. Consequently, to reconcile the need for authenticity with the operational requirements of the study, we developed a

Table 3B

Quality assurance of the VALGENT4 DNA samples: Amplifiable copy numbers at average and five different base pair length.

	Group	N ^a	Mean amplifiable copies	Ratio ^b	Standard deviation	95% Confidence Incidence		Min	Max	P-value ^c
						Lower	Upper			
Avg Copy	Valgent4	1,261	89,037	0.72	133,317	81,672	96,403	186	1,801,335	0.001
	Control2	182	123,942	1.00	114,870	107,141	140,743	3,000	879,343	
100 bp	Valgent4	1,261	34,782	1.06	84,948	30,089	39,475	369	1,633,069	0.757
	Control2	182	32,802	1.00	37,890	27,260	38,344	780	312,634	
200 bp	Valgent4	1,236	54,874	0.81	160,613	45,911	63,837	20	4,443,905	0.295
	Control2	181	67,766	1.00	104,680	52,413	83,119	569	795,228	
300 bp	Valgent4	1,261	156,201	0.77	380,175	135,198	177,205	193	8,254,550	0.124
	Control2	182	202,641	1.00	380,425	147,000	258,282	3,669	3,960,481	
400 bp	Valgent4	1,260	104,438	0.67	268,635	89,590	119,285	36	6,068,401	0.014
	Control2	182	156,585	1.00	268,209	117,356	195,813	3,726	2,904,381	
500 bp	Valgent4	1,249	93,544	0.59	396,191	71,550	115,537	48	8,947,189	0.031
	Control2	181	159,193	1.00	265,673	120,227	198,159	5,843	1,804,460	

^a Samples with insufficient number of SNPs were excluded, in addition a sample could have a valid SNPs number, but invalid result in one of the base pair groups, which mean that the numbers in the different base pair groups vary. In addition, samples with outlier values of above 10,000,000 were excluded from both sample sets.

^b The Control2 population was used as a reference.

^c The p-value is calculated using the one-way ANOVA test.

reconstitution protocol taking into account the changes introduced to the freshly collected VALGENT4 samples by routine cytology processing. The aim was to reestablish the residual cytology pelleted material to an approximate *per volume cellularity* corresponding to the original cellularity in un-manipulated clinical samples (Fig. 1). This reconstitution protocol was developed with the assistance of Becton, Dickinson and Company, the manufacturer of the SurePath medium, by reverse engineering the concentration and resuspension steps included in standard SurePath processing of original, fresh samples.

As the mean beta globin Ct value of the VALGENT4 samples (Ct = 24.8) was found to be 1.0 Ct above the Control-1 samples (Ct = 23.8), this suggests that the reconstituted panel samples had slightly lower cellularity than freshly collected, un-manipulated screening samples. However, in the context of the clinical positive-negative Onclarity HPV assay cut off (Ct 34.2), the reconstituted panel samples were fully acceptable as clinically, analytical material for HPV analysis. In combination with the ratio data from the ExomeQC analysis, this describes the resulting degree of dilution introduced to the reconstituted panel samples compared to the original sample cellularity. Moreover, the DNA fragmentation analysis showed that fresh, un-manipulated samples had more 400 and 500 bp DNA fragments compared to the reconstituted panel samples. This clear trend downwards between reconstituted samples and original samples in Control-2 was however not unexpected as all sorts of manipulation of cells and DNA leads to some level of deterioration. Whether this has implications for the analytical and clinical performance of HPV assays is unlikely, given that all included assays in the VALGENT4 protocol rely on amplification of less than 500bp HPV fragments (Table 1).

In conclusion, the two QA/QC methods employed indicate that the reconstituted samples used for the VALGENT evaluation contain slightly less analytical material compared to fresh, un-manipulated screening samples. However, the amount of analytical relevant material is reassuringly within the operational boundaries of all included VALGENT4 HPV assays (Table 1).

Besides the reconstitution protocol, VALGENT4 introduced central organization of DNA extraction for all assays running on this input material. This was chosen to eliminate any variation introduced at the level of DNA extraction methodology. Of the 11 assays evaluated in VALGENT4, only BD Onclarity, Seegene Anyplex HPV28/HPV II HR and Roche Cobas HPV assays required original SurePath LBC material as input material for analysis.

Until now, three VALGENT installments have been completed with several genotyping assays tested. The data from these VALGENT studies has showed that the VALGENT framework provides a good base for

cross-sectional clinical validation of HPV genotyping assays, by use of well annotated cervical screening samples collected in various screening programs around Europe [9]. A common feature of the VALGENT installments is the use of GP5 +/6 + PCR-EIA as the comparator assay to all other included assays. This brings the VALGENT in line with the stipulations within the international validation Guidelines of Meijer et al.

The first VALGENT study, VALGENT-1, included cervical cancer screening samples provided by the AML laboratory (Antwerp, Belgium) and included validation of four full genotyping assays [37–39]. The VALGENT-1 framework included SurePath samples from biobanked screening samples, and thus not freshly collected samples as included here. The second VALGENT study, VALGENT-2, used ThinPrep, involved the collection of samples from the Scottish Cervical Screening Programme via the Scottish HPV Reference laboratory and involved the evaluation of four assays with genotyping capability [40–43]. VALGENT-3 also included ThinPrep collected cervical cancer screening samples, this time from the Laboratory for Molecular Microbiology, Faculty of Medicine, University of Ljubljana, Slovenia [44] and evaluated 10 different HPV assays with different degrees of genotyping capability. VALGENT-1 [37–39] and VALGENT-2 [40–43] have both provided several peer-reviewed papers and publications from VALGENT-3 are underway [32,44]. Formal clinical validation of HPV assays for use in screening has primarily been undertaken on ThinPrep collected samples [24,26,28,30,41,45,46], with only one assay to date, the BD Onclarity, being validated on both ThinPrep [26] and SurePath collected samples [27]. Following this, it will be interesting to see how the HPV genotyping assays included in VALGENT4 perform on SurePath collected screening samples in this large-scale performance comparison.

Ethical and data inspection agency approvals

The study was approved by the Danish Data Inspection Agency J. No. AHH-2017-024, I-Suite: 05356. EU-GDPR compliant data handler agreement was established between the principal site Hvidovre Hospital and the Statistical analysis unit at Sciensano, Brussels. All collected samples were verified for non-compliance in the Danish human biological material in health research projects register (Væsanvendelsesregisteret).

Declaration of competing interest

VALGENT is an independent researcher induced research project,

set up by Sciensano, where manufacturers can have their HPV assays evaluated, under condition to provide equipment and kits and to cover costs for laboratory work and statistical analysis. Manufacturers cannot influence publication of manuscripts.

Jesper Bonde attended meetings with various HPV device manufacturers. JB has received honoraria from Hologic/Gen-probe, Roche, Qiagen, Genomica, and BD Diagnostics for lectures. Hvidovre Hospital has ongoing contracts with BD Diagnostics, Genomica, Self-Screen and EU Horizon2020.

Ditte Møller Ejegod attended meetings with various HPV device manufacturers.

Joakim Dillner has previously received research grants to his institution for research on HPV vaccines from an HPV vaccine manufacturer (Merck/SPMSD). The Karolinska University Hospital has ongoing contracts with several manufacturers of HPV diagnostics, including Roche and Genomica.

Daniëlle Heideman: minority stock portfolio in Self-screen B.V., a spin-off company of VU University Medical Center Amsterdam. Self-screen B.V. holds patents related to the work, and has developed and manufactured the HPV-Risk Assay. DAMH has been on the speaker's bureau of Qiagen and serves occasionally on the scientific advisory board of Pfizer and Bristol-Meyer Squibb.

Wim Quint: No conflicts of interest to declare.

Miguel Angel Pavon Ribas: received reagents for HPV testing at no cost from Roche, SeeGene, Qiagen and Genomica for research purposes.

Elizaveta Padalko: No conflicts of interest to declare.

Irene Kraus Christiansen: No conflicts of interest to declare.

Kate Cuschieri: No conflicts of interest to declare, KCs institution has received research funding and or associated gratis consumable from the following in the last 3 years: Hologic, Cepheid, Qiagen, Becton-Dickinson, Euroimmun, SelfScreen, LifeRiver, Genomica, GeneFirst.

Marc Arbyn and Lan Xu were supported by COHEAHR Network (grant No. 603019), coordinated by the Free University of Amsterdam (The Netherlands), funded by the 7th Framework Programme of DG Research and Innovation, European Commission (Brussels, Belgium).

Copenhagen University Hospital, Hvidovre: Reagents and instrumentation for testing Onclarity, CLART HPV4s and MassArray HPV test were received free of charge from the manufacturers for the duration of the testing. Limited co-funding for the project was obtained as part of the collaboration agreement with the manufacturers Becton, Dickinson and Company, Genomica SAU, and Agena Biotech GmbH.

Norwegian HPV Reference Lab: No external funding was obtained for this study.

Ghent University: Reagents for testing INNO-LIPA genotyping Extra II test was received free of charge from the manufacturer for the duration of the testing.

Infection and Cancer Laboratory, Barcelona: No external funding was obtained for this study

VU University Medical center: No external funding was obtained for this study.

Scottish HPV reference lab: Reagents and instrumentation for testing LifeRiver Venus, LifeRiver Harmonia, and GeneFirst HPV-MPA Genotyping Test were received free of charge from the manufacturers for the duration of the testing. Limited co-funding for the project was obtained as part of the collaboration agreement with the manufacturers Zhanghai Biotech (LifeRiver) and GeneFirst (HPV-MPA Genotyping Test).

DDL Diagnostic Laboratory: No external funding was obtained for this study.

Karolinska University hospital: No external funding was obtained for this study.

Sciensano, the institute where M.A. is employed, received support from Becton Dickinson, Genomica SAU, Agena Biotech GmbH, Zhanghai Biotech (LifeRiver), GeneFirst and FujiRebio for methodological and statistical work as described in the VALGENT Network (Arbyn et al., J Clin Virol 2016).

Author contribution

Design of Protocol: Jesper Bonde, Ditte Møller Ejegod, Lan Xu & Marc Arbyn.

Panel testing and reporting: Joakim Dillner, Daniëlle Heideman, Wim Quint, Miguel Angel Pavon Ribas, Elizaveta Padalko, Irene Kraus Christiansen, Kate Cuschieri, Ditte Møller Ejegod, Jesper Bonde.

Statistical analysis: Marc Arbyn & Lan Xu.

Writing of Manuscript: Jesper Bonde, Ditte Møller Ejegod, Marc Arbyn.

Editing of manuscript: All.

Decision to submit: All.

Acknowledgements

We wish to thank the staff of the Molecular Pathology Laboratory, Dept. Pathology, Hvidovre Hospital, Copenhagen, for their diligent work on assembling, aliquoting and processing the Valgent4 panel samples.

Moreover, we would like to thank all the participating company partners who engaged in this private-public research project.

References

- [1] V. Bouvard, R. Baan, K. Straif, Y. Grosse, B. Secretan, F. El Ghissassi, L. Benbrahim-Tallaa, N. Guha, C. Freeman, L. Galichet, et al., A review of human carcinogens—part B: biological agents, *Lancet Oncol.* 10 (4) (2009) 321–322.
- [2] Group IMw, IARCH monographs on the evaluation of carcinogenic risks of humans, Human Papillomavirus, IARC Press, Lyon, 2007, p. 90.
- [3] K.U. Petry, C. Liebrich, A. Luyten, M. Zander, T. Iftner, Surgical staging identified false HPV-negative cases in a large series of invasive cervical cancers, *Papillomavirus Res.* 4 (2017) 85–89.
- [4] H.A. Cubie, K. Cuschieri, Understanding HPV tests and their appropriate applications, *Cytopathology* 24 (5) (2013) 289–308.
- [5] M. Arbyn, G. Ronco, A. Anttila, C.J. Meijer, M. Poljak, G. Ogilvie, G. Koliopoulos, P. Naucler, R. Sankaranarayanan, J. Peto, Evidence regarding human papillomavirus testing in secondary prevention of cervical cancer, *Vaccine* 30 (Suppl. 5) (2012) F88–99.
- [6] M. Arbyn, J. Roelens, C. Simoons, F. Buntinx, E. Paraskevidis, P.P. Martin-Hirsch, W.J. Prendiville, Human papillomavirus testing versus repeat cytology for triage of minor cytological cervical lesions, *Cochrane Database Syst. Rev.* (3) (2013) CD008054.
- [7] G. Ronco, J. Dillner, K.M. Elfstrom, S. Tunesi, P.J. Snijders, M. Arbyn, H. Kitchener, N. Segnan, C. Gilham, P. Giorgi-Rossi, et al., Efficacy of HPV-based screening for prevention of invasive cervical cancer: follow-up of four European randomised controlled trials, *Lancet* 383 (9916) (2014) 524–532.
- [8] M. Poljak, J. Cuzick, B.J. Kocjan, T. Iftner, J. Dillner, M. Arbyn, Nucleic acid tests for the detection of alpha human papillomaviruses, *Vaccine* 30 (Suppl. 5) (2012) F100–106.
- [9] M. Arbyn, C. Depuydt, I. Benoy, J. Bogers, K. Cuschieri, M. Schmitt, M. Pawlita, D. Geraets, I. Heard, T. Gheit, et al., VALGENT: a protocol for clinical validation of human papillomavirus assays, *J. Clin. Virol.* 76 (Suppl. 1) (2016) S14–21.
- [10] M.J. Khan, P.E. Castle, A.T. Lorincz, S. Wacholder, M. Sherman, D.R. Scott, B.B. Rush, A.G. Glass, M. Schiffman, The elevated 10-year risk of cervical precancer and cancer in women with human papillomavirus (HPV) type 16 or 18 and the possible utility of type-specific HPV testing in clinical practice, *J. Natl. Cancer Inst.* 97 (14) (2005) 1072–1079.
- [11] X. Castellsague, K.A. Ault, F.X. Bosch, D. Brown, J. Cuzick, D.G. Ferris, E.A. Joura, S.M. Garland, A.R. Giuliano, M. Hernandez-Avila, et al., Human papillomavirus detection in cervical neoplasia attributed to 12 high-risk human papillomavirus genotypes by region, *Papillomavirus Res.* 2 (2016) 61–69.
- [12] J. Cuzick, C. Wheeler, Need for expanded HPV genotyping for cervical screening, *Papillomavirus Res.* 2 (2016) 112–115.
- [13] E.A. Joura, K.A. Ault, F.X. Bosch, D. Brown, J. Cuzick, D. Ferris, S.M. Garland, A.R. Giuliano, M. Hernandez-Avila, W. Huh, et al., Attribution of 12 high-risk human papillomavirus genotypes to infection and cervical disease, *Cancer Epidemiol. Biomarkers Prev.* 23 (10) (2014) 1997–2008.
- [14] V. Smelov, K.M. Elfstrom, A.L. Johansson, C. Eklund, P. Naucler, L. Arnheim-Dahlstrom, J. Dillner, Long-term HPV type-specific risks of high-grade cervical intraepithelial lesions: a 14-year follow-up of a randomized primary HPV screening trial, *Int. J. Cancer* 136 (5) (2015) 1171–1180.
- [15] S.K. Kjaer, C. Munk, J. Junge, T. Iftner, Carcinogenic HPV prevalence and age-specific type distribution in 40,382 women with normal cervical cytology, ASCUS/LSIL, HSIL, or cervical cancer: what is the potential for prevention? *Cancer Causes Control* 25 (2) (2014) 179–189.
- [16] M. Arbyn, M. Tommasino, C. Depuydt, J. Dillner, Are 20 human papillomavirus types causing cervical cancer? *J. Pathol.* 234 (4) (2014) 431–435.
- [17] D. Bzhalava, P. Guan, S. Franceschi, J. Dillner, G. Clifford, A systematic review of

- the prevalence of mucosal and cutaneous human papillomavirus types, *Virology* 445 (1–2) (2013) 224–231.
- [18] N. Li, S. Franceschi, R. Howell-Jones, P.J. Snijders, G.M. Clifford, Human papillomavirus type distribution in 30,848 invasive cervical cancers worldwide: variation by geographical region, histological type and year of publication, *Int. J. Cancer* 128 (4) (2011) 927–935.
 - [19] S.K. Kjaer, K. Frederiksen, C. Munk, T. Iftner, Long-term absolute risk of cervical intraepithelial neoplasia grade 3 or worse following human papillomavirus infection: role of persistence, *J. Natl. Cancer Inst.* 102 (19) (2010) 1478–1488.
 - [20] B. Serrano, S. de Sanjose, S. Tous, B. Quiros, N. Munoz, X. Bosch, L. Alemany, Human papillomavirus genotype attribution for HPV types 6, 11, 16, 18, 31, 33, 45, 52 and 58 in female anogenital lesions, *Eur. J. Cancer* 51 (13) (2015) 1732–1741.
 - [21] K. Kavanagh, K.G. Pollock, K. Cuschieri, T. Palmer, R.L. Cameron, C. Watt, R. Bhatia, C. Moore, H. Cubie, M. Cruickshank, et al., Changes in the prevalence of human papillomavirus following a national bivalent human papillomavirus vaccination programme in Scotland: a 7-year cross-sectional study, *Lancet Infect. Dis.* 17 (12) (2017) 1293–1302.
 - [22] W.K. Huh, K.A. Ault, D. Chelmow, D.D. Davey, R.A. Goulart, F.A. Garcia, W.K. Kinney, L.S. Massad, E.J. Mayeaux, D. Saslow, et al., Use of primary high-risk human papillomavirus testing for cervical cancer screening: interim clinical guidance, *Obstet. Gynecol.* 125 (2) (2015) 330–337.
 - [23] M. Arbyn, L. Xu, F. Verdoodt, J. Cuzick, A. Szarewski, J.L. Belinson, N. Wentzensen, J.C. Gage, M.J. Khan, Genotyping for human papillomavirus types 16 and 18 in women with minor cervical lesions: a systematic review and meta-analysis, *Ann. Intern. Med.* 166 (2) (2017) 118–127.
 - [24] C.J. Meijer, J. Berkhof, P.E. Castle, A.T. Hesselink, E.L. Franco, G. Ronco, M. Arbyn, F.X. Bosch, J. Cuzick, J. Dillner, et al., Guidelines for human papillomavirus DNA test requirements for primary cervical cancer screening in women 30 years and older, *Int. J. Cancer Suppl.* 124 (3) (2009) 516–520.
 - [25] M. Arbyn, P.J. Snijders, C.J. Meijer, J. Berkhof, K. Cuschieri, B.J. Kocjan, M. Poljak, Which high-risk HPV assays fulfil criteria for use in primary cervical cancer screening? *Clin. Microbiol. Infect.* 21 (9) (2015) 817–826.
 - [26] D. Ejegod, I. Serrano, K.S. Cuschieri, W.A. Nussbaumer, L.M. Vaughan, A.S. Ahmad, J. Cuzick, J. Bonde, Clinical validation of the BD onclarity HPV assay using a non-inferiority test, *Med. Microbiol. Diagn.* (2013) S3, <https://doi.org/10.4172/2161-0703.S3-003>.
 - [27] D. Ejegod, F. Bottari, H. Pedersen, M.T. Sandri, J. Bonde, The BD onclarity HPV assay on samples collected in SurePath medium meets the international guidelines for human papillomavirus test requirements for cervical screening, *J. Clin. Microbiol.* 54 (9) (2016) 2267–2272.
 - [28] D.A. Heideman, A.T. Hesselink, J. Berkhof, F. van Kemenade, W.J. Melchers, N.F. Daalmeijer, M. Verkuijten, C.J. Meijer, P.J. Snijders, Clinical validation of the cobas 4800 HPV test for cervical screening purposes, *J. Clin. Microbiol.* 49 (11) (2011) 3983–3985.
 - [29] A.T. Hesselink, D.A. Heideman, J. Berkhof, F. Topal, R.P. Pol, C.J. Meijer, P.J. Snijders, Comparison of the clinical performance of PapilloCheck human papillomavirus detection with that of the GP5+/6+-PCR-enzyme immunoassay in population-based cervical screening, *J. Clin. Microbiol.* 48 (3) (2010) 797–801.
 - [30] A.T. Hesselink, C.J. Meijer, M. Poljak, J. Berkhof, F.J. van Kemenade, M.L. van der Salm, M. Bogaarts, P.J. Snijders, D.A. Heideman, Clinical validation of the Abbott RealTime High Risk HPV assay according to the guidelines for human papillomavirus DNA test requirements for cervical screening, *J. Clin. Microbiol.* 51 (7) (2013) 2409–2410.
 - [31] A.T. Hesselink, R. Sahli, J. Berkhof, P.J. Snijders, M.L. van der Salm, D. Agard, M.C. Bleeker, D.A. Heideman, Clinical validation of Anyplex II HPV HR Detection according to the guidelines for HPV test requirements for cervical cancer screening, *J. Clin. Virol.* 76 (2016) 36–39.
 - [32] L. Xu, A. Ostrbenk, M. Poljak, M. Arbyn, Assessment of the Roche linear array HPV genotyping test within the VALGENT framework, *J. Clin. Virol.* 98 (2018) 37–42.
 - [33] C. Eklund, O. Forslund, K.L. Wallin, J. Dillner, Continuing global improvement in human papillomavirus DNA genotyping services: The 2013 and 2014 HPV LabNet international proficiency studies, *J. Clin. Virol.* 101 (2018) 74–85.
 - [34] L.S.A. Muhr, C. Eklund, J. Dillner, Towards quality and order in human papillomavirus research, *Virology* 519 (2018) 74–76.
 - [35] K. Pedersen, S. Fogelberg, L.H. Thamsborg, M. Clements, M. Nygard, I.S. Kristiansen, E. Lynge, P. Sparen, J.J. Kim, E.A. Burger, An overview of cervical cancer epidemiology and prevention in Scandinavia, *Acta Obstet. Gynecol. Scand.* 97 (7) (2018) 795–807.
 - [36] D.M. Ejegod, J. Junge, M. Franzmann, B. Kirschner, F. Bottari, M. Sideri, M.T. Sandri, J. Bonde, Clinical and analytical performance of the BD Onclarity HPV assay for detection of CIN2+ lesions on SurePath samples, *Papillomavirus Res.* 2 (2016) 31–37.
 - [37] M. Schmitt, C. Depuydt, I. Benoy, J. Bogers, J. Antoine, M. Pawlita, M. Arbyn, group Vs, Viral load of high-risk human papillomaviruses as reliable clinical predictor for the presence of cervical lesions, *Cancer Epidemiol. Biomarkers Prev.* 22 (3) (2013) 406–414.
 - [38] M. Schmitt, C. Depuydt, I. Benoy, J. Bogers, J. Antoine, M. Arbyn, M. Pawlita, Group VS, Prevalence and viral load of 51 genital human papillomavirus types and three subtypes, *Int. J. Cancer* 132 (10) (2013) 2395–2403.
 - [39] M. Schmitt, C. Depuydt, I. Benoy, J. Bogers, J. Antoine, M. Arbyn, M. Pawlita, V.S. Group, Multiple human papillomavirus infections with high viral loads are associated with cervical lesions but do not differentiate grades of cervical abnormalities, *J. Clin. Microbiol.* 51 (5) (2013) 1458–1464.
 - [40] K. Cuschieri, D. Geraets, J. Cuzick, L. Cadman, C. Moore, D. Vanden Broeck, E. Padalko, W. Quint, M. Arbyn, Performance of a cartridge-based assay for detection of clinically significant human papillomavirus (HPV) infection: lessons from VALGENT (validation of HPV genotyping tests), *J. Clin. Microbiol.* 54 (9) (2016) 2337–2342.
 - [41] K. Cuschieri, D.T. Geraets, C. Moore, W. Quint, E. Duvall, M. Arbyn, Clinical and analytical performance of the onclarity HPV assay using the VALGENT framework, *J. Clin. Microbiol.* 53 (10) (2015) 3272–3279.
 - [42] D.T. Geraets, K. Cuschieri, M.N. de Koning, L.J. van Doorn, P.J. Snijders, C.J. Meijer, W.G. Quint, M. Arbyn, Clinical evaluation of a GP5+/6+-based luminex assay having full high-risk human papillomavirus genotyping capability and an internal control, *J. Clin. Microbiol.* 52 (11) (2014) 3996–4002.
 - [43] I. Heard, K. Cuschieri, D.T. Geraets, W. Quint, M. Arbyn, Clinical and analytical performance of the PapilloCheck HPV-Screening assay using the VALGENT framework, *J. Clin. Virol.* 81 (2016) 6–11.
 - [44] N.J. Polman, A. Ostrbenk, L. Xu, P.J.F. Snijders, C. Meijer, M. Poljak, D.A.M. Heideman, M. Arbyn, Evaluation of the clinical performance of the HPV-Risk assay using the VALGENT-3 panel, *J. Clin. Microbiol.* 55 (12) (2017) 3544–3551.
 - [45] A. Boers, R. Wang, L. Slagter-Menkema, B.M. van Hemel, H. Ghyssaert, A.G. van der Zee, G.B. Wisman, E. Schuur, Clinical validation of the Cervista HPV HR test according to the international guidelines for human papillomavirus test requirements for cervical cancer screening, *J. Clin. Microbiol.* 52 (12) (2014) 4391–4393.
 - [46] D.A. Heideman, A.T. Hesselink, F.J. van Kemenade, T. Iftner, J. Berkhof, F. Topal, D. Agard, C.J. Meijer, P.J. Snijders, The Aptima HPV assay fulfills the cross-sectional clinical and reproducibility criteria of international guidelines for human papillomavirus test requirements for cervical screening, *J. Clin. Microbiol.* 51 (11) (2013) 3653–3657.